



PERIODONTAL

Sickle Scalers	D5
Curettes/Universal	D13
Curettes/Area-Specific	D16
Hoe, File & Chisel Scalers	D26
Sidekick® Sharpener & Sharpening Stones	D29
Educational Aids	D31
Miscellaneous	D32



PERIODONTAL SET-UPS

LARGE PROPHY SET-UP

DESCRIPTION	SUGGESTION	PAGE
Cassette	IM416x*	A3
Mirror/Handle	MIR5/3, MH6	B15
Explorer	EXD11/126	B4
Probe	PCP12/11.5	B8
Sickle Scalars	SH6/79, S204S9	D5, D7
Universal Curette	SC13/149	D14
Area Specific Curettes	SG1/29, SG7/89, SG11/129, SG13/149	D16, D17
Implant Scaler	IMPLACARE	D25
Dressing Pliers	DP2	B13
Sharpening Stones	SS4, SS299	D29, D30
Pumice Cup	PCUPM	D32

Other Recommended Products:

Prophy Angle, Suction Tip, Aspirating Syringe, A/W Syringe Tip, IMS Cassette Clip, and Color-Coded Rings.

ROOT PLANING SET-UP

DESCRIPTION	SUGGESTION	PAGE
Cassette	IM416x*	A3
Mirror/Handle	MIR5/3, MH6	B15
Explorers	EXD11/12A6, EXTU17/236	B3, B4
Probe	PCP15/11.5	B8
Sickle Scalars	SH6/76, S204SD6	D5, D7
Root Planing Curettes	SAS1/29, SAS7/89, SRPG11/129, SRPG13/149, SRPG15/169, SG17/189	D17, D20, D21
Universal Curette	SM17/189	D11
Dressing Pliers	DP2	B13
Sharpening Stones	SS3C, SS299	D29, D30
Pumice Cup	PCUPM	D32

Other Recommended Products:

Prophy Angle, Suction Tip, Aspirating Syringe, A/W Syringe Tip, IMS Cassette Clip, and Color-Coded Rings.

SMALL PROPHY SET-UP

DESCRIPTION	SUGGESTION	PAGE
Cassette	IM508x*	A2
Mirror/Handle	MIR4/3, MH6	B15
Explorer	EXD2A	B3
Probe	PCP12/11.5	B8
Sickle Scalars	SH6/79, S204S9	D5, D7
Universal Curettes	SC13/149, SBH1/29, SBH5/69	D14
Pumice Cup	PCUPM	D32

Other Recommended Products:

Prophy Angle, Suction Tip, Color-Coded Rings, Sharpening Stone SS3C.

PEDODONTIC/SEALANT/IMPLANT SET-UP

DESCRIPTION	SUGGESTION	PAGE
Cassette	IM508x*	A2
Mirror/Handle	MIR4/3, MH6	B15
Expro	XP23/11.56	B11
Dressing Pliers	DP18	B13
Sickle Scalars	SH6/79, S204S9	D5, D7
Curette	SC13/149	D14
Implant Scaler	IMPLACARE	D25
Pumice Cup	PCUPM	D32

Other Recommended Products:

Prophy Angle, Suction Tip, Color-Coded Rings, Sharpening Stone SS3C.

See Surgical Set-Ups for Periodontal Surgery Set-Up, GP Periodontal Surgery Set-Up, and Suture Removal Set-Up.

x* See Section A for exact product code number





HANDLE OPTIONS

Most popular instruments available in all handle options.

DOUBLE-END

#9
EverEdge®
Technology



#8
ResinEight™



#7
Satin Steel®
Colours™



#6
Satin Steel®



#4
Round



#2
Octagonal



SINGLE-END

#8
ResinEight™



#7
Satin Steel®
Colours™



#6
Satin Steel®



#24
Round



#1
Octagonal



HOW TO USE THIS SECTION

11/12
Scaler

Instrument Name & Pattern

SM11/127

Pictured Part Code

#2, #4, #6, #7,
#8, #9

Available handle designs
for this pattern.

See index for all available part codes of a specific pattern.

DID YOU KNOW?

Carpal Tunnel Syndrome Prevention

Neurologists recommend alternating instrument handle sizes as one means of reducing stress. Larger diameter handles (#6, #7, #8) help lighten instrument grasp. Using a combination of various handle sizes plus a more relaxed grasp can help lessen the severity of the symptoms of Carpal Tunnel Syndrome. *

* Source: Gerwatowski, L.J., McFall, D.B., Stach, D.: Carpal Tunnel Syndrome; Risk Factors and Preventive Strategies for the Dental Hygienist. Journal of Dental Hygiene, February 1992.



EVEREDGE® TECHNOLOGY

STAYS SHARPER LONGER.

Revolutionary new stainless steel alloy is super-durable to stay sharper longer. That means less time sharpening and less hand fatigue.

PROPRIETARY PROCESS

Proprietary, state-of-the-art technology in metallurgy, heat treatment and cryogenic processing ensures that the superior edge retention and wear characteristics of EverEdge Technology will last the entire life of the instrument.

The improved sharpness of EverEdge Technology is not a superficial coating – the long-lasting wear is present throughout the entire instrument tip. EverEdge Technology scalars can be sharpened again and again for your best instrument value.

Throughout my career, I was always the delegated "instrument sharpener". I accepted this role eagerly because I knew how important it was for the patient and for the clinician to use a sharp edge during instrumentation. Life is too short and your clinical day is filled with important obligations, the Everedge is a lifesaver.

Cheryl Thomas, RDH

Hu-Friedy EverEdge instruments made a positive change in my practice of Dental Hygiene. The variety of EverEdge instruments available allows access to any surface. A consistently sharp blade reduces hand fatigue and increases patient comfort. My hand instrument set-ups now consist entirely of Hu-Friedy EverEdge instruments!

Lory Laughter, RDH, BS

FINELY-HONED FINISH

EverEdge Technology Scalars and Curettes feature a finely-honed finish to the instrument blade for super sharp edges and efficient scaling.

DIAMOND KNURL

Unique diamond knurl pattern with large diameter handle for a confident grasp. Optimal weight for reduced hand fatigue.

COMFORT ZONE™ CONE WITH SIGNATURE SERIES™ GRIP

The cushioned grip provides a smooth transition and increased comfort. Distinctive Ocean Blue grip for easy identification.





EVEREDGE® SCALERS & PERIODONTAL CURETTES HANDLE GUIDE

GRACEY STANDARD

SG1/29
SG3/49
SG5/69
SG7/89
SG9/109
SG11/129
SG13/149
SG11/149
SG12/139
SG15/169
SG17/189

GRACEY RIGID

SG1/2R9
SG3/4R9
SG5/6R9
SG7/8R9
SG11/12R9
SG13/14R9
SG15/16R9
SG17/18R9

GRACEY AFTER FIVE®

SRPG1/29
SRPG3/49
SRPG5/69
SRPG7/89
SRPG11/129
SRPG13/149
SRPG15/169

GRACEY AFTER FIVE RIGID

SRP1/2R9
SRP3/4R9
SRP5/6R9
SRP7/8R9
SRP11/12R9
SRP13/14R9
SRP15/16R9

GRACEY MINI FIVE®

SAS1/29
SAS3/49
SAS5/69
SAS7/89
SAS11/129
SAS13/149
SAS15/169

GRACEY MINI FIVE RIGID

SAS1/2R9
SAS3/4R9
SAS5/6R9
SAS7/8R9
SAS11/12R9
SAS13/14R9
SAS15/16R9

GRACEY MICRO MINI FIVE

SMS1/29
SMS7/89
SMS11/129
SMS13/149

PINK RIBBON CURETTE

SYME9

SICKLE SCALERS

S204S9
S204SD9
SH5/339
SH6/79
SIUFW2049
SJ30/339
SJ31/329
SJ34/359
SM0/009
SN1359
SN1379
SU15/309
SU15/339

NEVI® SERIES

SCNEVI19
SCNEVI29
SCNEVI39
SCNEVI49
SNEV1/H59

UNIVERSAL CURETTES

BARNHART

SBH1/29
SBH1/2R9
SBH5/69
SBH5/6R9

COLUMBIA

SC2R/2L9
SC4R/4L9
SC13/149

YOUNGER-GOOD

SYG7/89

McCALL

SM13/149
SM13/14S9
SM17/189

LANGER

SL1/29
SL3/49
SL5/69
SL17/189

UC-RADCLIFF

SUCR3/49

SIGNATURE SERIES™ GRIPS

Gray	SSG1
White	SSG2
Orange	SSG3
Brown	SSG4
Yellow	SSG5
Red	SSG6
Green	SSG7
Blue	SSG8
Black	SSG9
Pink	SSG10
Purple	SSG11
Ocean Blue*	SSG12
16 Assorted	SSG0
80 Assorted	SSGMULTI

All single color packages contain 16 pieces.

The assorted color packages contain your choice of 16 or 80 pieces.

**Default color for all EverEdge items.*



SCALERS & PERIODONTAL CURETTES HANDLE GUIDE

	SINGLE-END HANDLES					DOUBLE-END HANDLES					
	1	24	6	7	8	2	4	6	7	8	9
Gracey						.	.	.	.	.	.
After Five Gracey						.	.	.	.	.	.
Mini Five Gracey						.	.	.	.	.	.
Ball		.					.				
Barnhart						.	.	.	.	.	.
Bunting						.	.	.			
Bunting File	.						.	.			
Cattoni		.	.				.	.	.	.	
Columbia						.	.	.	.	.	.
Goldman-Fox						.	.	.			
Hirschfeld File						.	.	.			
Indiana University	.				.	.	.	.	.		
Jacquette	.	.	.			.	.	.	.	.	
Kramer-Nevins						.	.				
Langer						.	.	.	.	.	
McCall		.				.	.	.	.	.	.
Morse		.				.	.	.	.	.	.
Odont	.					.	.	.			
Nevi								.	.	.	.
Ratcliff							.		.		
Quétin							.	.			
SH6/7, SH5/33, S204, S204S, S204SD						.	.	.	.	.	.
SN135, SN137, SN137M							.	.	.	.	.
Taylor		.				.	.	.			
U15		.	.	.	.	.	.	.	.	.	.
U15/30, U15/33							.	.	.	.	.
Vision Curvettes							.	.			
Younger-Good	.					.	.	.	.	.	.

If the instrument tip/handle combination you desire is not shown, we will be glad to provide you with a quote for its preparation.
 Note: Instrument tips and handles are not sold separately.



SICKLE SCALERS

Designed for the removal of supragingival calculus.



H6/H7
Hygienist

Super-thin, contra-angle design. Very effective for interproximal areas. Shank angulation provides accessibility to all anterior and premolar surfaces.

SH6/79

#2, #4, #6, #7,
#8, #9




6/7
Goldman

Shortened blade version of H6/H7

SHG6/76

#2, #4, #6,
#7, #8

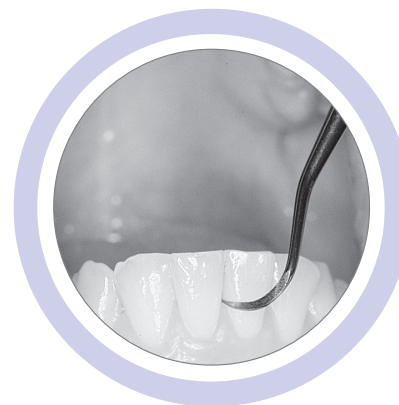

SH6/7 on Anterior Teeth



U15
Towner

SU156


#24, #6, #7, #8



H5/U15
Hygienist-
Towner

SH5/U156


#2, #4, #6

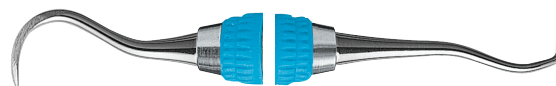

U15/30
Towner-
Jacquette

SU15/309

#2, #4, #6, #7,
#8, #9

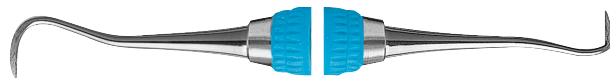


H5/33
Hygienist-
Jacquette

SH5/339

#2, #4, #6, #7,
#8, #9

Super-thin sickle scaler combined with a straight sickle blade. Very effective for scaling interproximal and cervical areas of anterior teeth.

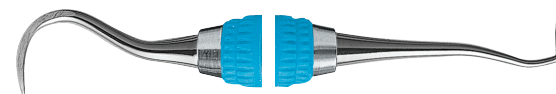

U15/33
Towner-
Jacquette

SU15/339

#2, #4, #6, #7,
#8, #9


DID YOU KNOW?

Some scalers and curettes are available in an American Dental Classic style that is slightly more rigid in design, with a $\frac{5}{32}$ " working end. They are offered exclusively in ResinEight™ handles. Contact your Hu-Friedy representative for specific part code availability.



SICKLE SCALERS

NEVI® SCALERS

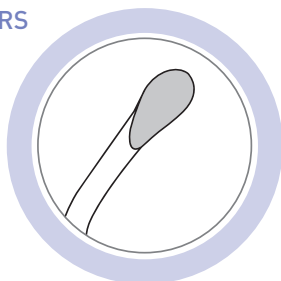

Close-up of Nevi 1 disk end



Nevi 1
Anterior

An elongated disk end for anterior lingual stain and calculus removal paired with a super thin, anterior curved sickle.

SCNEVI16

#7, #8, #9



Nevi 1/H5
Anterior

An elongated disk end for anterior lingual stain and calculus removal paired with a larger H5 sickle – a thin curved blade for interproximal and cervical anterior areas.

SNEV1/H59

#6, #7, #8, #9



Nevi 2
Posterior

A super thin posterior curved sickle scaler. A unique contra-angle design great for interproximal reach and ergonomic hand positioning.

SCNEVI29

#6, #7, #8, #9



Nevi 3
Posterior

Ultra slender, dual cutting blades provide effective placement in tight, narrow contacts and interproximal areas and precise access for effective removal of light-to-moderate deposits.

SCNEVI39

#9




Nevi 4
Posterior


Tapered dual cutting blades provide precise access for effective deposit removal. Increased shank rigidity provides for effective removal of burnished, heavy and tenacious deposits

SCNEVI49

#9

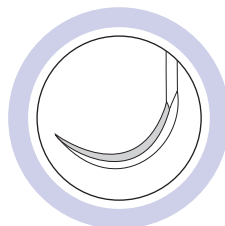

THE CURVED SICKLE SCALER
Two cutting edges on a curved blade that end in a sharp point

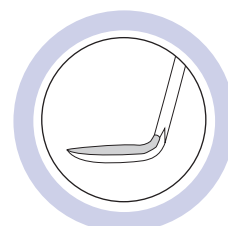

THE STRAIGHT SICKLE SCALER
Two cutting edges on a straight blade that end in a sharp point. Also known as a Jacquette Scaler

JACQUETTE SCALERS*



30/33
Jacquette

Combination standard and miniature anterior Jacquette

SJ30/339

#2, #4, #6, #7,
#8, #9



31/32
Jacquette

Posterior standard Jacquette

SJ31/329

#2, #4, #6, #7,
#8, #9



34/35
Jacquette

Posterior miniature Jacquette

SJ34/359

#2, #4, #6, #7,
#8, #9



H3/H4
Jacquette

SH3/46

#2, #4, #6



14/15
Posterior
Jacquette

SJ14/156

#2, #4, #6

DID YOU KNOW?

* Jacquette #30 - same as 1Y
Jacquette #31 - same as 2Y
Jacquette #32 - same as 3Y

Miniature Jacquette #33 same as 1YS
Miniature Jacquette #34 same as 2YS
Miniature Jacquette #35 same as 3YS



SICKLE SCALERS



204

S2046

#2, #4, #6



204
IUFW

SIUFW2049

#2, #4, #6, #7,
#8, #9



204S

For posterior supra-gingival interproximal scaling.

S204S9

#2, #4, #6, #7,
#8, #9



204SD

Contra angle design offers greater access to distal surfaces.

S204SD9

#2, #4, #6, #7,
#8, #9



1
Springfield
Anterior

SSP16

#24, #6



2/3
Springfield

SSP2/36

#2, #4, #6



SICKLE SCALERS

POSTERIOR



135

SN1359

#4, #6, #7,
#8, #9



REM 3/4
Remington

SRM3/46

#2, #4, #6



129

S1296

#2, #4, #6



T2/T3
Taylor

ST2/36

#2, #4, #6,
#7, #8

Available single-ended (ST2 and ST26).



SICKLE SCALERS

ANTERIOR



137

Small anterior sickle paired with an anterior universal curette.

SN1379

#4, #6, #7, #8, #9



137
Mini Five®

SN137M6

#6, #7, #8



0/00
Morse

Miniature blade for scaling crowded or overlapping anterior teeth. Excellent for stain removal.

SM0/009

#2, #4, #6, #7, #8, #9



OD1/33
Odont/
Jacquette

SOD1/336

#2, #4, #6



128
Nebraska

SN128

#24, #6, #7, #8



USC-128

SUSC128

#24, #6



WH-2
Whiteside

SWH2

#24



YG-15
Younger-
Good

SYG156

#1, #6



SICKLE SCALERS

SICKLE SCALERS



CI-1

SCI1

#24, #6




152
Cumine

SCM152

#2

Use spoon end for stain removal.
Use sickle end for supragingival scaling.

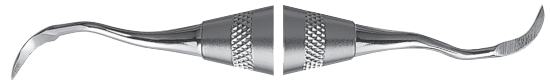

CI-2/3

SCI2/36

#2, #4, #6, #7

BALL SCALERS



B-1

SB1


#24



B-4

SB4

#24

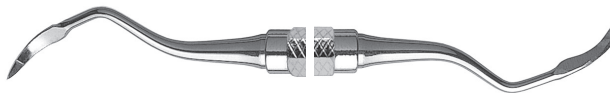

B-2/3

SB2/3

#2




B-5/6

SB5/6

#2, #4

CRANE-KAPLAN

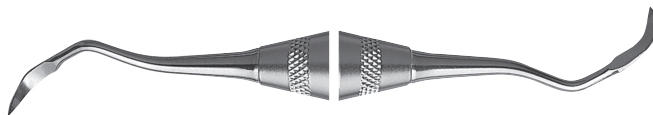

CK-6

SCK66

#2, #4, #6



CK-6D

SCK6D

#2, #4

Angle is not as acute and the blade is curved in a different plane than the SCK6. Overall blade size is smaller.



SCALER/CURETTE* SETS


McCALL SET



11/12
Scaler

SM11/127


#2, #6, #7, #8

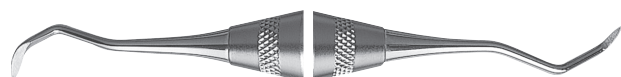

11A/12A
Scaler

SM11/12A6

#2, #4, #6

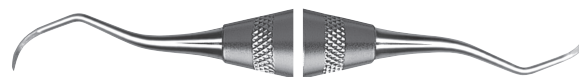

11/12S
Modified
McCall

S11/12S6

#4, #6




13S/14S
Pointed
Curette

SM13/14S9

#2, #4, #6, #7,
#8, #9

Round tip available (SM13/14).

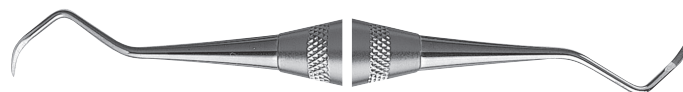

MTC 13/14
Pointed
Curette

SMT13/146

#4, #6

Longer blade (1mm) than M13/14S.



17S/18S
Curette

SM17/189

#2, #4, #6, #7,
#8, #9

SM17/18S is the same as SM17/18.
Pointed tip available (SM17/18PT).



19/20
Curette

SM19/20

#2

* Note: This publication spells curette as per Merriam-Webster's Collegiate Dictionary, Eleventh Edition, ©2005 Merriam-Webster, Inc.



SCALER/CURETTE SETS

GOLDMAN-FOX SET




GF-1
Universal
Scaler

SGF16

#4, #6

For removing heavy supragingival and subgingival calculus. Useful on proximal surfaces and interproximal areas during periodontal surgical procedures.

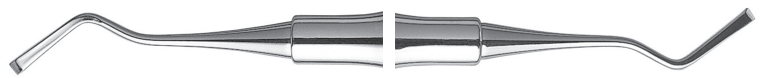

GF-5 Hoe
Buccal/
Lingual

SGF5

#4

For removing heavy calculus deposits.




GF-2
Curette

SGF26

#4, #6

Contour design for removal of calculus on anterior teeth and premolars.

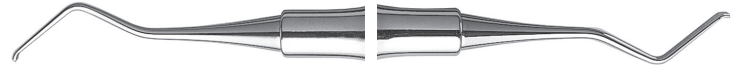

GF-6
Hoe

SGF6

#4

For mesial/distal surfaces of the posterior teeth or facial/lingual surfaces of the anterior teeth.



GF-3
Curette

SGF36

#4, #6

Same blade surface as GF-2, specifically designed for the proximal surface of the posterior teeth.



GF-21

SGF21

#4, #6

Sharp spoon on one end for debridement in curettage procedures. Small scaler on the other end can be used with either a circumferential or vertical stroke.




GF-4
Curette

SGF46

#4, #6

Additional bend in the shank for easier calculus removal in the molar region.

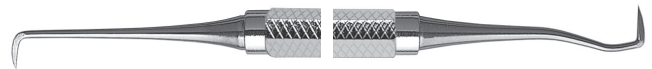

GF-21S/33
Combination

SGF21S/33

#2, #4

Small end Goldman-Fox scaler with a #33 straight miniature Jacquette scaler.

CATTONI SET



103/106

S103/1066

#4, #6



107/108

S107/1086


#2, #4, #6, #7, #8



104/105
Hoe/Chisel

S104/1056

#4, #6

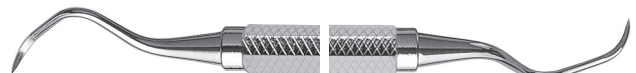

107/108S
Modified

S107/108S

#2, #4

Shorter shank than the S107/108.



CURETTES/UNIVERSAL

LANGER CURETTES



1/2 Langer (shown)	1/2 Langer Rigid	1/2 Langer After Five®	1/2 Langer Mini Five®
SL1/29	SL1/2R6	SL1/2AF6	SL1/2MF6
#4, #6, #7, #8, #9	#4, #6, #7, #8	#4, #6	#4, #6

Langer 1/2 for mandibular posterior.



3/4 Langer (shown)	3/4 Langer Rigid	3/4 Langer After Five®	3/4 Langer Mini Five®
SL3/49	SL3/4R6	SL3/4AF6	SL3/4MF6
#4, #6, #7, #8, #9	#4, #6, #7, #8	#4, #6	#4, #6

Langer 3/4 for maxillary posterior.



5/6 Langer (shown)	5/6 Langer Rigid	5/6 Langer After Five®	5/6 Langer Mini Five®
SL5/69	SL5/6R6	SL5/6AF6	SL5/6MF6
#4, #6, #7, #8, #9	#4, #6	#4, #6	#4, #6



Langer 5/6 for maxillary and mandibular anterior.



17/18 Langer (shown)	17/18R Langer Rigid
SL17/189	SL17/18R6
#4, #6, #7, #8, #9	#4, #6

Langer 17/18 for maxillary and mandibular posterior.

LANGER APPLICATION CHART

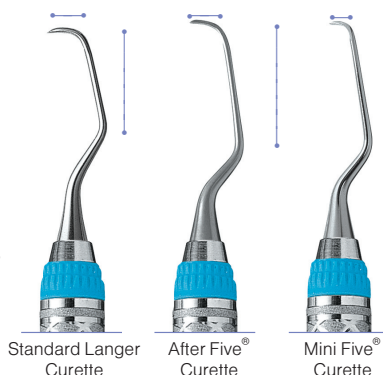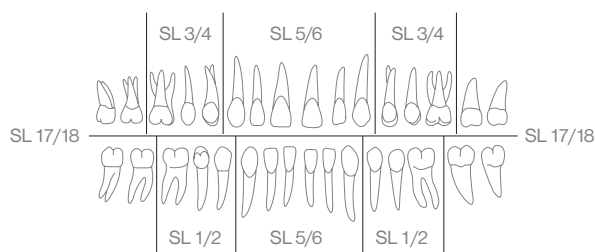

- After Five and Mini Five terminal shanks are elongated by 3mm for access into deep periodontal pockets and root surfaces of 5mm or more.
- Blade is 10% thinner than standard Langer to ease gingival insertion and reduce tissue distention.
- Mini Five has a reduced blade, half the length of the After Five or standard Langer, for better adaptation in narrow pockets and furcations.


DID YOU KNOW?

Langer curettes are designed with the shank design of a Gracey curette combined with a universal blade. These instruments can scale all surfaces.




CURETTES/UNIVERSAL

Blade designed with two cutting edges and a rounded toe for use on all surfaces.

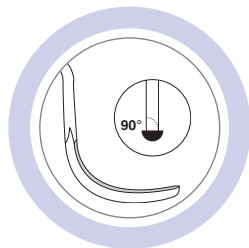

Blade angulation of a
Universal Curette

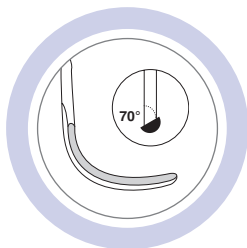

Blade angulation of a
Gracey Curette



IU13/14
Indiana
University

SIU13/146

#2, #4, #6



IU17/18
Indiana
University

SIU17/186

#2, #4, #6, #7, #8




1/2
Bunting

Curette-Chisel combination

SBN1/2

#2, #4, #6

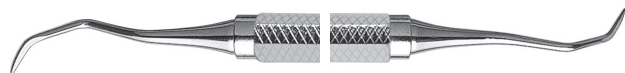

5/6
Bunting

SBN5/6

#2



YG7/8
Younger-Good

SYG7/89

#2, #4, #6, #7,
#8, #9

BARNHART



1/2
Barnhart
(shown)

SBH1/29


#2, #4, #6, #7,
#8, #9



1/2
Barnhart
Rigid

SBH1/2R8

#2, #4, #6, #7,
#8, #9

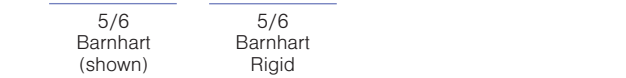

5/6
Barnhart
(shown)

SBH5/69

#2, #4, #6, #7,
#8, #9



5/6
Barnhart
Rigid

SBH5/6R8

#2, #4, #6, #7,
#8, #9

COLUMBIA



2R/2L
Columbia

Anterior application

SC2R/2L9

#2, #4, #6, #7,
#8, #9



4R/4L
Columbia

Posterior application

SC4R/4L9


#2, #4, #6, #7,
#8, #9



13/14
Columbia
(shown)

SC13/149

#2, #4, #6, #7,
#8, #9

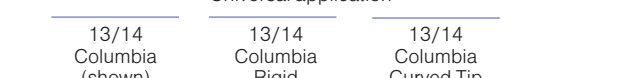

13/14
Columbia
Rigid

SC13/14R8

#2, #4, #6,
#7, #8




13/14
Columbia
Curved Tip

SC13/14C

#2, #4



CURETTES/UNIVERSAL

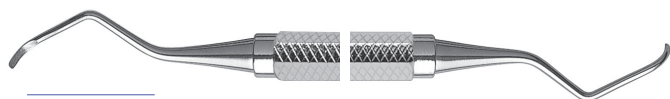

13S/14S
Bates

SBT13/14S

#2



HU1/2
Hutchinson

SHU1/2

#2



3/4
Rule [23-24]

SR3/4


#2, #4, #6



UC3/4
Ratcliff

SUCR3/49

#6, #8, #9

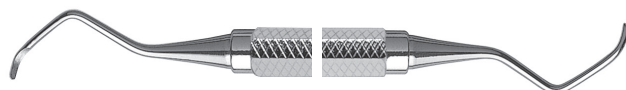

411/412

S411/412

#2, #4



10/11
Loma Linda
(shown)

SLL10/116

#6, #7, #8

10/11
Loma Linda
Rigid


SLL10/11R6

#6, #8

10/11
Loma Linda
Slim

SLL10/11S6


#6, #8

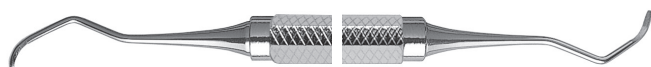

413/414

S413/414


#2, #4

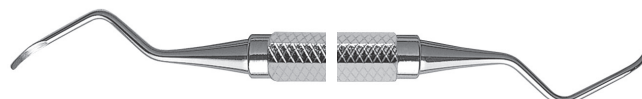

1/2
Mallery

SMAL1/2

#2, #4

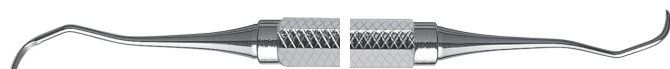


31/32
Goldman

SHG31/32

#2

Curette version of 31/32 Jacquette.

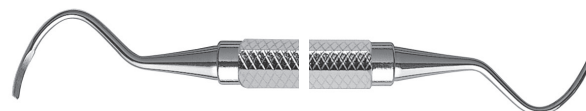

ME-1
[10-11]

SME1


#2




GRACEY STANDARD CURETTES

Designed to adapt to a specific area or tooth surface.


- Area-specific to allow for deep scaling, root planing and periodontal debridement.
- Offset blade provides a perfect working angulation for the tooth surface.
- Offset blade angulation means only the lower cutting edge of each blade is used.

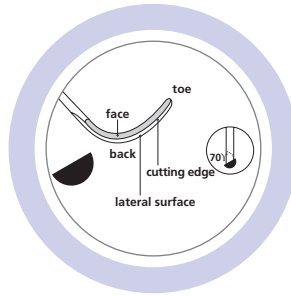

Blade angulation of a Gracey Curette

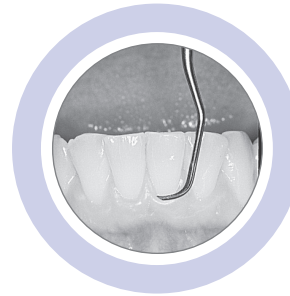

Anterior SG 1/2

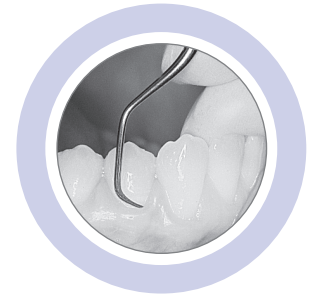

Posterior SG 7/8



1/2
(shown)

1/2
Rigid

Slight contra-angle for anterior incisors and canines.

SG 1/29

#2, #4, #6, #7,
#8, #9

SG 1/2R9

#2, #4, #6, #7,
#8, #9



3/4
(shown)

3/4
Rigid

Short contra-angle for anterior incisors and canines.

SG 3/49

#2, #4, #6, #7,
#8, #9

SG 3/4R9

#2, #4, #6, #7,
#8, #9



5/6
(shown)

5/6
Rigid

5/6
Extra Rigid

Medium contra-angle for anteriors and premolars.

SG 5/69

#2, #4, #6, #7,
#8, #9

SG 5/6R9

#2, #4, #6, #7,
#8, #9

SG 5/6ER6

#4, #6, #7, #8



7/8
(shown)

7/8
Rigid

7/8
Extra Rigid

Medium contra-angle for premolars and molars, facial and lingual surfaces.

SG 7/89

#2, #4, #6, #7,
#8, #9

SG 7/8R9

#2, #4, #6, #7,
#8, #9

SG 7/8ER6


#4, #6, #7, #8

DID YOU KNOW?

Extra Rigid Gracey curettes are ideal for tenacious calculus removal, gross scaling, and initial debridement. Blades, however, are identical in size to standard Gracey curettes.




GRACEY CURETTES

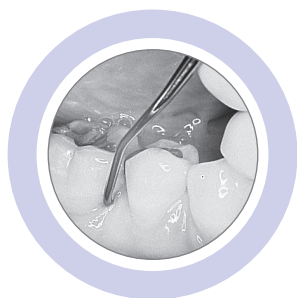

Posterior SG11/12



Posterior SG13/14

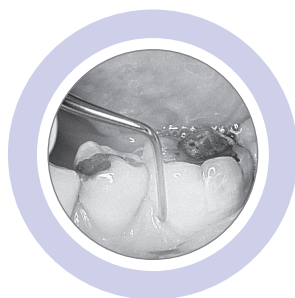

Posterior SG15/16



Posterior SG17/18



9/10⁺
(shown)

SG9/109
#2, #4, #6, #7,
#8, #9

9/10
Rigid

SG9/10R9
#2, #4, #6,
#7, #8

Long contra-angle for molar buccal and lingual surfaces. More pronounced angle than 7/8 for less accessible root surfaces.



11/12⁺
(shown)

SG11/129
#2, #4, #6, #7,
#8, #9

11/12
Rigid

SG11/12R9
#2, #4, #6, #7,
#8, #9

11/12
Extra Rigid

SG11/12ER8
#4, #6, #7, #8

Angulated to reach mesial surfaces of posterior teeth.



13/14⁺
(shown)

SG13/149
#2, #4, #6, #7,
#8, #9

13/14
Rigid

SG13/14R9
#2, #4, #6, #7,
#8, #9

13/14
Extra Rigid

SG13/14ER8
#4, #6, #7, #8

Angulated to reach distal surfaces of posterior teeth.



15/16⁺
(shown)

SG15/169
#2, #4, #6, #7,
#8, #9

15/16
Rigid

SG15/16R9
#2, #4, #6, #7,
#8, #9

Same shank bends as 13/14, but blade positioned to reach mesial posterior surfaces.



17/18

SG17/189
#2, #4, #6, #7,
#8, #9

17/18
Rigid

SG17/18R9
#2, #4, #6, #7,
#8, #9

Accentuated angles improve access to distal posterior surfaces. Also has smaller blade and slightly longer terminal shank.

DID YOU KNOW?

The Rigid Gracey has a thicker, stronger, and less flexible shank than the standard Gracey. This allows for moderate to heavy calculus removal without having to utilize additional heavy scalers such as sickles and hoes.



GRACEY CURETTES

Designed to adapt to a specific area or tooth surface.

GRACEY MESIAL-DISTAL CURETTES

These combinations of working ends are designed to complete an entire facial or lingual sextant without changing instruments, increasing efficiency.



11/14
Mesial/Distal
(shown)

SG11/146

#2, #4, #6,
#7, #8

11/14
Mesial/Distal
Rigid

SGR11/14R

#2, #4



12/13
Mesial/Distal
(shown)

SG12/136

#2, #4, #6,
#7, #8

12/13
Mesial/Distal
Rigid

SGR12/13R

#2, #4

GRACEY PROPHY CURETTES

The blades are the same size as the Gracey Standard Curettes, but the shank is shorter and more rigid for removal of supra-gingival deposits.



P1/2

SGP1/2

#2, #7



P3/4

SGP3/4

#2, #7



P5/6

SGP5/6

#2, #7



P7/8

SGP7/8

#2, #7



THE CLAYTON GRACEY STORY

In the late 1940s, an innovative dentist from the University of Michigan had a vision. Clayton Gracey wanted to "give every dentist the possibility to treat even the deepest and least accessible periodontal pockets simply and without traumatic stretching of the gingiva."

Dr. Gracey brought his idea to Hugo Friedman, the founder of Hu-Friedy, who supported his concept. Together they developed a series of 14 single-ended area-specific instruments for the successful removal of deposits from root surfaces.

Dr. Gracey's designs laid the groundwork for the new instruments of today. Currently, Gracey curettes are available as double-ended instruments with multiple handle choices.


The Gracey curettes have come a long way since the original design. Hu-Friedy currently manufactures a variety of styles and patterns, such as After Five® curettes with a longer terminal shank and thinned blade, Mini Five® curettes with After Five features and a shorter blade, and the 15/16 and 17/18 patterns. An optional Rigid shank is available on all of the Gracey designs.

Hu-Friedy is continually looking for new ways to improve our existing products. Our success is dedicated to the pioneers in dentistry, including Dr. Clayton Gracey, who helped to make a vision become a reality and a part of every dental office today.



GRACEY CURETTES

GRACEY CURETTES – DESIGN COMPARISONS

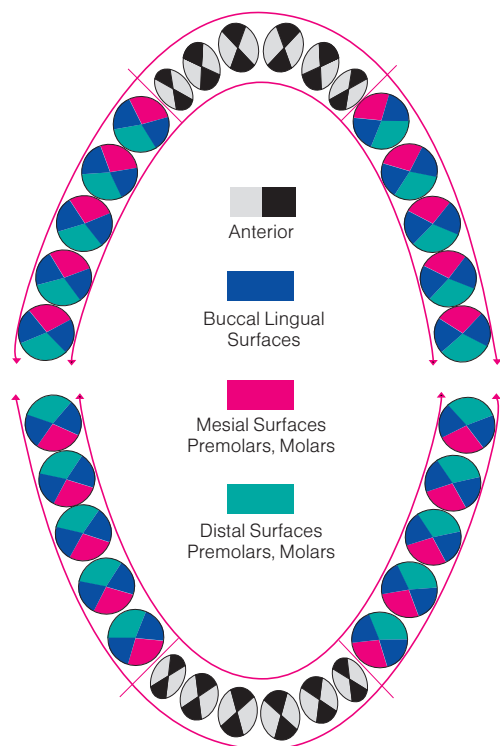

Refer to the chart to match instruments with area of application by using the color-coded diagram above.

-  1/2, 3/4, 5/6
-  7/8, 9/10
-  11/12, 15/16
-  13/14, 17/18*

	Shank Design & Diameter	Blade Length	Blade Width	Available Patterns & Area of Use see chart to match color
Standard (Finishing)	Standard	Standard	Standard	1/2, 3/4, 5/6 7/8, 9/10 11/12, 15/16, 13/14, 17/18*
Rigid	Standard design, increased shank diameter	Standard	Standard	1/2, 3/4, 5/6 7/8, 9/10 11/12, 15/16, 13/14, 17/18*
After Five	Longer terminal shank, standard diameter	Standard	Decreased by 10%	1/2, 3/4, 5/6 7/8 11/12, 15/16, 13/14
Rigid After Five	Longer terminal shank, increased diameter	Standard	Decreased by 10%	1/2, 3/4, 5/6 7/8 11/12, 15/16, 13/14
Mini Five	Longer terminal shank, standard diameter	Decreased by 50%	Decreased by 10%	1/2, 3/4, 5/6 7/8 11/12, 15/16, 13/14
Rigid Mini Five	Longer terminal shank, increased diameter	Decreased by 50%	Decreased by 10%	1/2, 3/4, 5/6 7/8 11/12, 15/16, 13/14
Micro Mini Five	Longer terminal shank, standard diameter	Decreased by 50%	Decreased by 20% compared to Mini Five	1/2, 7/8 11/12, 13/14

DID YOU KNOW?


Hu-Friedy Thinning Service
Most scalers and curettes (except After Five® and Mini Five®) can be thinned to create an extra fine blade. Add the following to the original part code:
TSDE to double-end or TSSE to single-end

* The 17/18 is a unique pattern, having a longer terminal shank and slightly shorter blade.



GRACEY CURETTES – AFTER FIVE®

Designed to adapt to a specific area or tooth surface.

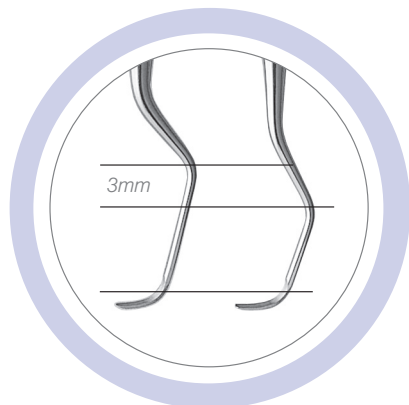

Terminal shank redesigned to be 3mm longer than a standard Gracey curette



Gracey After Five curette reaching deep into a pocket

AFTER FIVE® CURETTES

- Terminal shank elongated by 3mm for access into deep periodontal pockets and root surfaces of 5mm or more and better clearance around the crown.
- Blade thinned by 10% to ease gingival insertion and reduce tissue distention.
- Area-specific design patterned after the standard Graceys offers exact adaptation and control.
- Available in Rigid and Extra Rigid designs for greater strength when needed.

1/2	3/4	5/6	7/8	11/12	13/14	15/16	11/14 Mesial/Distal	12/13 Mesial/Distal
SRPG1/29	SRPG3/49	SRPG5/69	SRPG7/89	SRPG11/129	SRPG13/149	SRPG15/169	SRPG11/149	SRPG12/139
#2, #4, #6, #7, #8, #9* on all After Five patterns.								
Rigid	Rigid	Rigid	Rigid	Rigid	Rigid	Rigid	Rigid	Rigid
SRP1/2R9	SRP3/4R9	SRP5/6R9	SRP7/8R9	SRP11/12R9	SRP13/14R9	SRP15/16R9	SRP11/14R4	SRP12/13R4
#2, #4, #6, #7, #8, #9* on all After Five Rigid patterns except where noted.							#4	#4
		Extra Rigid	Extra Rigid	Extra Rigid	Extra Rigid			
		SRP5/6E6	SRP7/8E6	SRP11/12E6	SRP13/14E6			
#6, #7, #8 on all After Five Extra Rigid patterns.								




GRACEY CURETTES – MINI FIVE®

MINI FIVE® CURETTES

- Terminal shank elongated by 3mm for access into deep periodontal pockets and root surfaces of 5mm or more.
- Reduced blade, half the length of the After Five® or standard Gracey, for better adaptation in narrow pockets and furcations.
- Blade is 10% thinner than standard Gracey to ease gingival insertion and reduce tissue distention.
- Area-specific design patterned after the standard Graceys offers exact adaptation and control.
- Available in Rigid designs for greater strength when needed.



Mini Five reaching into a deep pocket on a narrow root

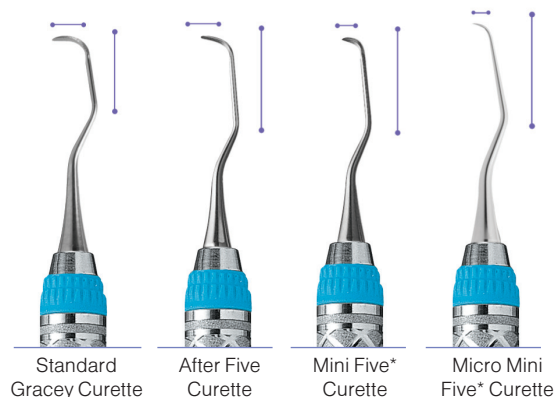


Standard Gracey Curette

After Five Curette

Mini Five* Curette

Micro Mini Five* Curette

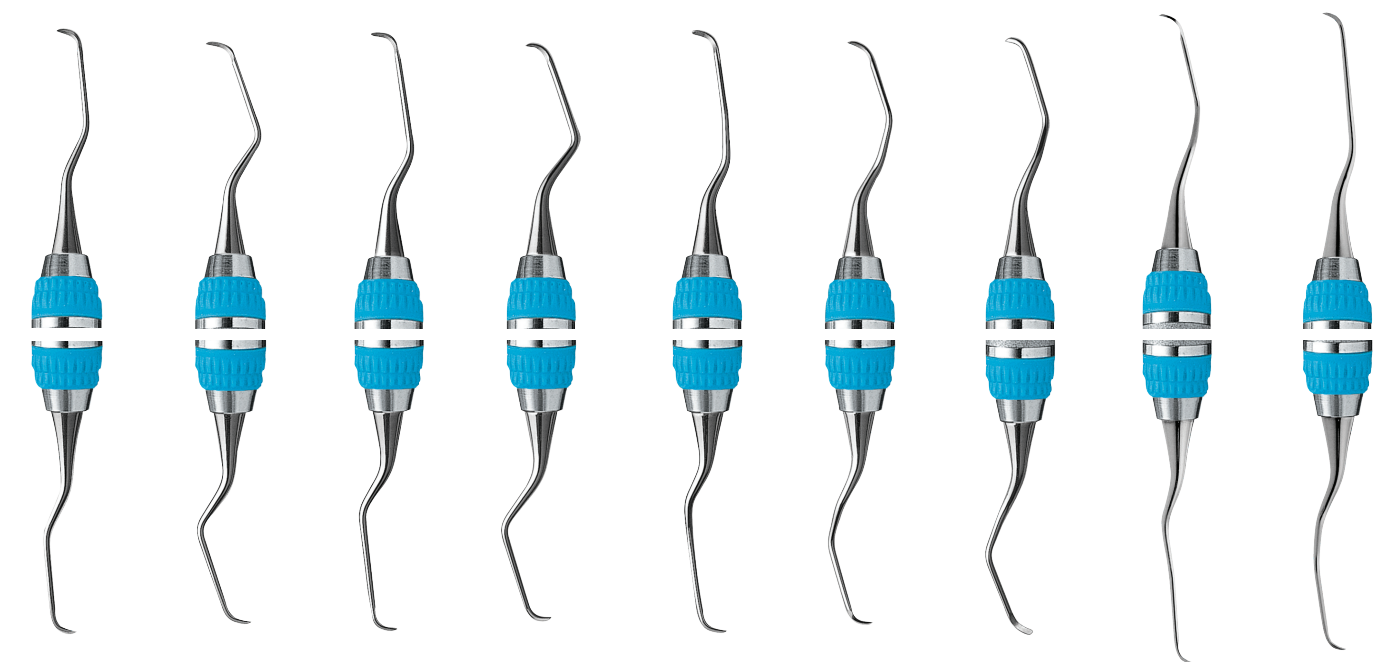

1/2

3/4

5/6

7/8

11/12

13/14

15/16

11/14
Mesial/Distal

12/13
Mesial/Distal

SAS1/29

SAS3/49

SAS5/69

SAS7/89

SAS11/129

SAS13/149

SAS15/169

SAS11/149

SAS12/139

#2, #4, #6, #7, #8, #9 on all Mini Five patterns.

Rigid

Rigid

Rigid

Rigid

Rigid

Rigid

Rigid

SAS1/2R9

SAS3/4R9

SAS5/6R9

SAS7/8R9

SAS11/12R9

SAS13/14R9


SAS15/16R9

#2, #4, #6, #7, #8, #9 on all Mini Five Rigid patterns.



GRACEY CURETTES – MICRO MINI FIVE®

- Elongated terminal shank for access into deep periodontal pockets and root surfaces of 5mm and more.
- Reduced blade, half the length of the After Five or standard Gracey, for better adaptation in narrow pockets and furcations.
- Blade is 20% thinner than a Mini Five Gracey to further reduce tissue dissention and ease gingival insertion.
- Area-specific design patterned after the standard Gracey curettes offers exact adaptation and control.
- Revolutionary EverEdge® Technology stays sharper longer for less time sharpening and reduced hand fatigue.

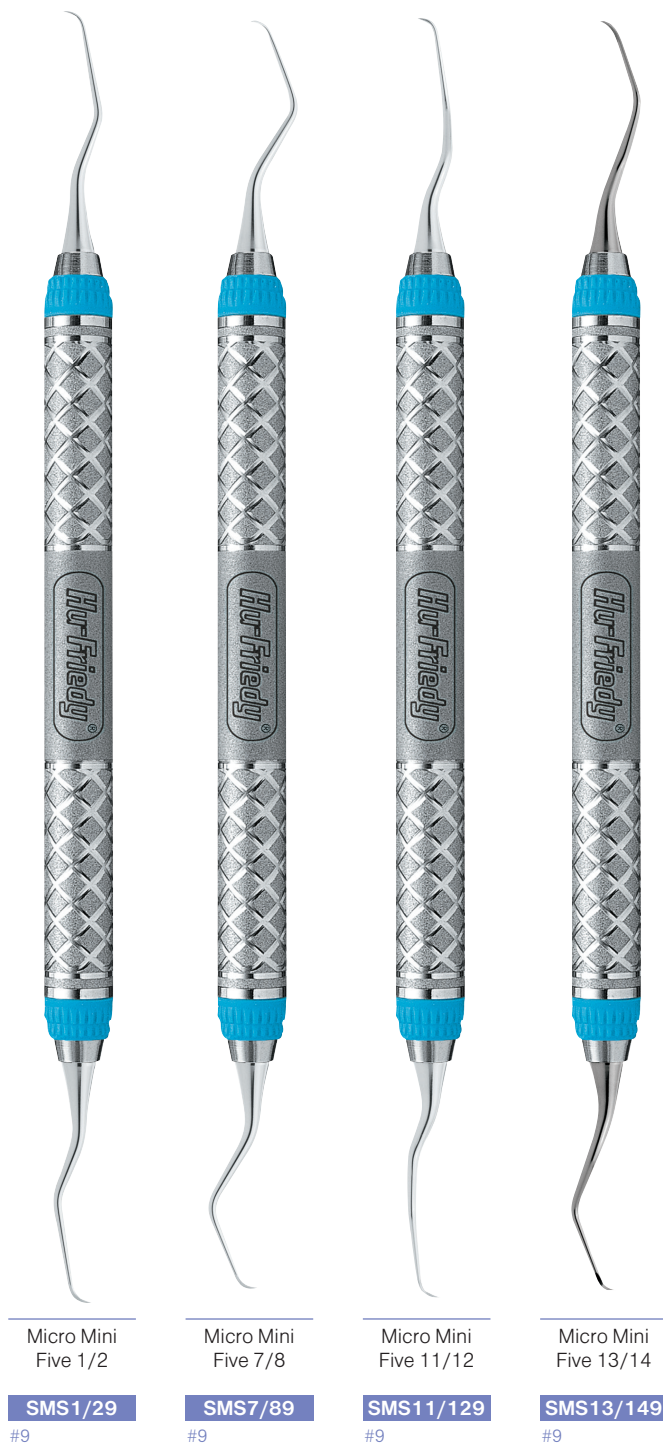




CURETTES/AREA-SPECIFIC

Designed to adapt to a specific area or tooth surface.

VISION CURVETTES®

- 50% shorter, curved blades.
- 5mm and 10mm shank markings.
- Blade identification mark (+) on handle.
- Identification bands on handle. One band on 11/12, two bands on 13/14 (on #4 handle only).

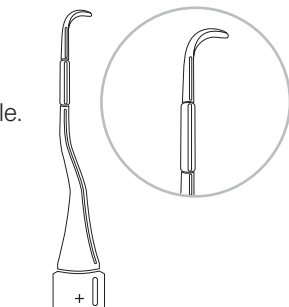


Curvette
11/12

For mesial surfaces of molars.
Improved adaptation into furcations.

SGC11/12

#4, #6

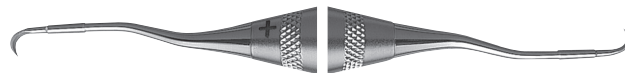


Curvette
Sub-Zero

Long shank reaches deep into the
sulcus, adapting to the facial and
lingual root surfaces.

SGCSUB-06

#4, #6

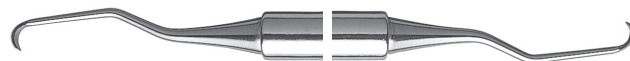

Curvette
1/2

For anterior and premolar surfaces,
and interproximal areas.

SGC1/2

#4, #6



Curvette
13/14


For distal surfaces of molars.
Improved adaptation into furcations.

SGC13/14

#4, #6

TURGEON MODIFIED GRACEY CURETTES

- Blade is designed to produce a different cross section of blade angulation.
- Modified shape provides a sharper blade.
- Modified shape is easy to sharpen.
- Blades are narrow for ease of insertion.

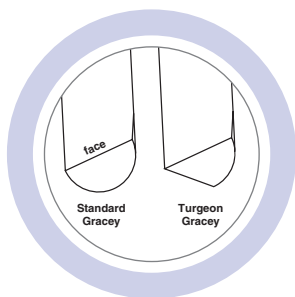

*Angulation of a Standard Gracey
vs. a Turgeon Gracey*



5/6
Modified
Gracey

STG5/66

#2, #6



11/12
Modified
Gracey

STG11/126

#2, #6



13/14
Modified
Gracey

STG13/146

#2, #6



13/14S
Modified
McCall

STM13/14S6

#2, #6



17/18S
Modified
McCall

STM17/18S6

#2, #6




CURETTES/AREA-SPECIFIC


DEBRIDEMENT CURETTES*

These unique curettes are used after ultrasonic scaling for gentle removal of residual deposits and for smoothing the root surface.


IDEAL FOR FURCATIONS, DEVELOPMENTAL GROOVES AND LINE ANGLES.

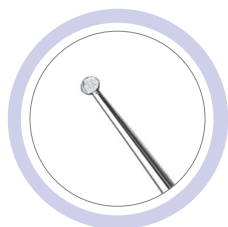

Entire edge of the blade is a cutting edge enabling a push or pull stroke in all directions – vertical, horizontal or oblique.

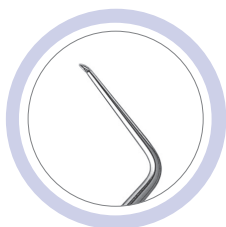


Shanks are extended for easy access into deep pockets.

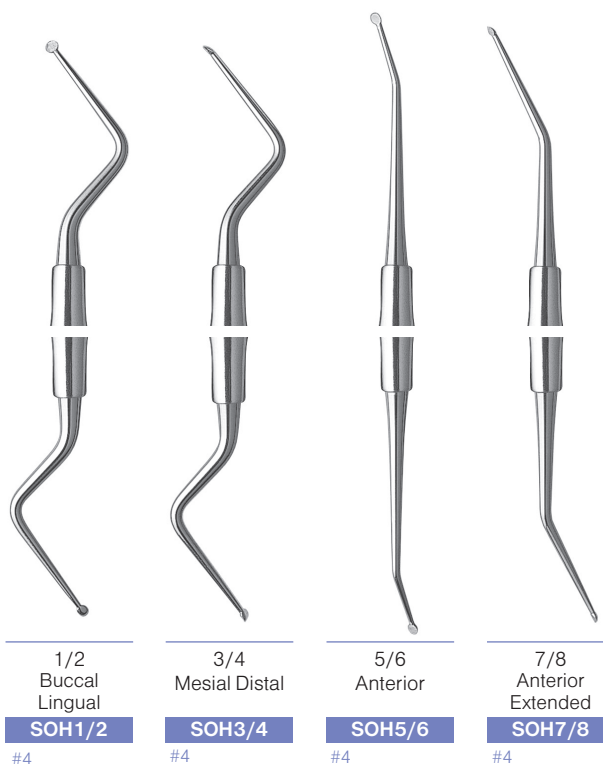

FURCATION CURETTES

Designed for scaling root concavities and furcations.

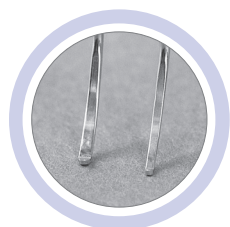

Comparison of 1.3mm and 0.9mm curette blades



DID YOU KNOW?

The DeMarco Curettes are excellent for stain removal on anterior lingual surfaces.




IMPLACARE™


For implant maintenance. Used to remove plaque and calculus without damaging titanium abutments.




Implacare 1 Curette Kit cleaning an abutment.
This curette utilizes the Universal Columbia 4R/4L tip design.

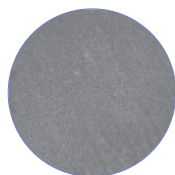

High power magnification of a titanium abutment.

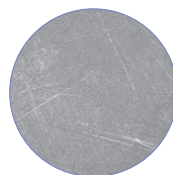

After 25 strokes with a plastic material that contains abrasives.

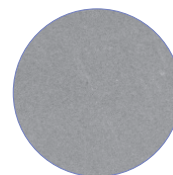

After 25 strokes with Plasteel™ Implacare tips.

IMPLANT MAINTENANCE KIT

Implant Maintenance Care Starter Kit

IMPKIT

#6

- 1 Implacare scaler handle
- 1 set of 204S scaler tips
- 1 set of Columbia 4R/4L tips
- 1 set of H6/H7 Scaler tips
- 1 Colorvue probe handle
- 1 set of PCP12 tips (3-6-9-12)
- 1 set of PCP11 tips (3-6-8-11)
- 1 set of Williams tips (1-2-3-5-7-8-9-10)
- 1 set of UNC12 (1-2-3-4-5-6-7-8-9-10-11-12)
- Includes evidence-based implant maintenance protocol



Implant Maintenance Kit 1

IMPLACAR16

#6
Columbia 4R/4L Plasteel tips, and Handle.
(10 pairs)

Implant Maintenance Kit 2

IMPLACAR26

#6
204S Plasteel tips, and Handle.
(10 pairs)

Implant Maintenance Kit 3

IMPLACAR36

#6
H6/H7 Plasteel tips, and Handle.
(10 pairs)

Implant Maintenance Kit A

IMPLACARA6

#6
Assorted tips, and Handle.
(12 pairs)

Columbia 4R/4L Plasteel Tips

IC4R/4L

(10 pairs)

204S Plasteel Tips

I204S

(10 pairs)

H6/H7 Plasteel Tips

IH6/7

(10 pairs)

Assorted Tips

IMPA

(12 pairs)

Implacare Satin Steel Handle

IMPHDL6

#6

REFILL PACKAGES


DID YOU KNOW?

The Implacare Curette Tips are made of Plasteel and are disposable. Plasteel is a high grade resin that does not scratch titanium and should be autoclaved or chemiclaved before use.



HOE SCALERS

Used to dislodge heavy supramarginal calculus deposits.
Use a pull stroke.

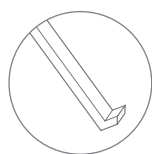

XS-15
Towner

SXS156

#24, #6

For anterior
linguals



6/7
Orban
Mesial/Distal
SO6/7

#2, #6



8/9
Orban Buccal/
Lingual
SO8/9

#2, #6



U15/13K

SU15/13K

#4



WV1/2

SWV1/2

#2




7K

S7K

#24

FILE SCALERS


Used to crush and remove heavy calculus deposits.
Use a pull stroke.

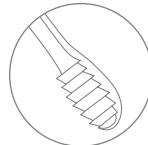

3/7
Hirschfeld
Buccal/Lingual

FH3/76

#2, #4, #6

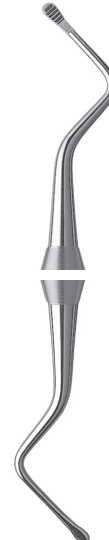

10/11
Orban Buccal/
Lingual
FO10/11

#2, #4, #6



Orban
Mesial/Distal
FO12/13

#2, #4, #6



5/11
Hirschfeld
Mesial/Distal

FH5/11

#2, #4, #6




9/10
Hirschfeld
Anterior
FH9/10

#2, #4, #6



CHISEL SCALERS

Used to dislodge bridges of calculus on anterior mandibular teeth. Use a push stroke.

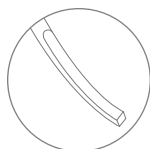

62
Cushing

SCU62

#1



61
Cushing

SCU61

#1



Zerfing

SZ6

#24, #6

DARBY-PERRY CHISELS



DP5K/6K

SDP5/6K

#2



DP9K/11K

SDP9/11K

#2




DIAMOND COATED FILE SCALERS

Used for definitive calculus removal. Use multidirectional strokes.



DIAMONDTEC™ FILE SCALERS

- The perfect size and shape to enhance definitive scaling in furcations and root depressions.
- Double-ended designs allow for more efficient scaling.
- Large diameter handle provides an ergonomic grip and helps reduce the symptoms of CTS*.
- Signature Series grips make instruments easy to identify.

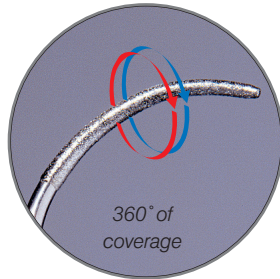

SDCN7



SDCM/D7

KRAMER-NEVINS PERIODONTAL SURGICAL KIT

KRAMER-NEVINS CURETTES

The Kramer-Nevins modified Gracey curette series is designed for root planing and subgingival curettage. Features slightly modified blade curvatures, thinner profiles, and extra length to help gain access to difficult areas without losing curette "bite." Most commonly employed during initial pre-surgical preparation, and with all maintenance patients. Also extremely useful during surgery.



#5/6
Kramer-Nevins
Curette

SKN5/6



#7/8
Kramer-Nevins
Curette

SKN7/8



#11/12
Kramer-Nevins
Curette

SKN11/12



#13/14
Kramer-Nevins
Curette

SKN13/14



SIDEKICK® SHARPENER

The Hu-Friedy Sidekick sharpener is designed for routine maintenance sharpening of scalers and curettes. Correct routine sharpening helps extend instrument life.



FEATURES

- The Sidekick sharpener is cordless, for use anywhere in the dental office.
- Instrument Guide Channels and Vertical Backstops are "template-like" guides that help control blade angulation, providing consistent results.
- The Sidekick sharpener is compact and requires minimal space.

Sidekick Sharpening Kit

SDKKIT

Kit Includes: Sidekick Unit, User Manual & CD-ROM, Ceramic Sharpening Stone, Plastic Test Stick, Magnifying Lens and 2 AA Batteries.

REPLACEMENT PARTS & ACCESSORIES

Guideplate Screw	SDKGPSCREW	Stone Clip	SDKCLIP
Guideplate	SDKGPLATE	Plastic Test Sticks	PTS
Ceramic Stone	SDKSTONEC	Magnifying Lens	LENS

SHARPENING STONES



No.1
I-Stone
Fine Grit °

SSF1

4" x 1" x 1/2"
(10 x 2.5 x 1.3cm)



No.3
Ceramic
Medium Grit °

SS3C

3" x 1 1/4" x 1/4"
(7.6 x 3.2 x .6cm)



No.4
Arkansas Flat
Hard
Fine Grit °

SS4

4 1/8" x 1" x 3/8"
(10 x 2.5 x .6cm)



No.5
Arkansas Bates
Hard Fine Grit °

SS5

4" x 1 1/2" x 1/2"
(10 x 3.8 x 1.3cm)



No.6
I-Stone
Wedge Shape
Medium Grit °

SS6

4 1/2" x 1 7/8" x 3/8"
down to 1/8" (11.4
x 4.8 x .9 down to
.3cm)



No.6A
Arkansas
Thick Wedge Hard,
Fine Grit °

SS6A

4" x 1 7/8" x 3/8" down
to 1/8" (11.4 x 4.8 x
.9 down to .3cm)



AS-22
Arkansas
Super Hard
Super Fine Grit °

SSAS22

2 1/4" x 7/8" x 3/16" (5.7
x 2.2 x .5cm)



MS-22
I-Stone
Medium Grit °

SSMS22

2 1/4" x 7/8" x 3/16" (5.7
x 2.2 x .5cm)



Dual Grit
Stone

SSDGRIT

Features a coarse
ceramic side for
recontouring and a light
ceramic side
for finishing and honing.
4" x 5/8" x 1".

SHARPENING TIP

For information on "How to Sharpen"
refer to pages K3-K6 in the catalog
reference section.

W - Lubricate stone with water
O - Lubricate stone with oil
D - No lubrication needed

Sidekick US Patent No.6,254,464



SHARPENING STONES

Stones illustrated 50% size.



299
Arkansas Conical
Super-Fine Grit[®]

SS299

3½" x 5/16" (8.9 x .8cm) base
Hard, for sharpening
facial surface of
curette.



No.2
Arkansas
Cylindrical
Super-Fine Grit[®]

SS2

4" x ¼" (10 x .6cm)
For sharpening
facial surface of
curette.



Ceramic
Cylindrical
Fine Grit[®]

SSKC

4" x 5/16" (10 x .31cm)



309
I-Stone
Fine Grit[®]

SS309

4" x ½" to ¼" x ¼" to
⅛" (10 x 1.3 x .6 x
.6 x .3cm)



File Sharpener
Tanged
Carbon Steel

FST

4¼" (10.8cm) long



Hu-Friedy File
Sharpener
Carbon Steel

FS

5½" (14cm) long



Plastic Test Stick
Autoclavable

PTS

3" x ¼" (7.6 x .6cm)
Available in
packages of 6.



Luxating Elevator
Sharpening Stone
Fine Grit[®]

SSLUX

For use with luxator
elevators.

COMPOSITION STONES FOR THE HANDPIECE

Coarse. Made of white aluminum oxide Shown full size.



#301^w

SS301

¼" x ¾" (6 x 10mm)
diameter



#303^w

SS303

¼" x ½" (6 x 13mm)
diameter



#304^w

SS304

¼" x ¾" (6 x 16mm)
diameter



#305^w

6 x 19mm

SS305

¼" x ¾" (6 x 19mm)
diameter



#300^w

6 x 6mm

SS300

¾" x ¼" (6 x 6mm)
at base



#302^w

3 x 3mm

SS302

¼" x ¼" (3 x 3mm)
diameter



Sharpen-Ez™
Sharpening Oil
1-oz. Bottle
29.5ml

SSO

Pre-sharpening
application to stone
floats away metal
particles and prevents
damaging build-up.
Also decreases
friction heat.



Arkansas
Stone Kit

SSKITA

- #4 Flat Arkansas Stone
- #299 Conical Arkansas Stone
- Plastic Test Sticks (5)
- Sharpen-EZ Stone Oil
- Magnifying Lens



Ceramic
Stone Kit

SSKITC

- #3 Ceramic Stone
- Cylindrical Ceramic Stone
- Plastic Test Sticks (5)
- Magnifying Lens



Sharpening
Magnifying
Lens

LENS


W - Lubricate stone with water
O - Lubricate stone with oil
D - No lubrication needed



EDUCATIONAL AIDS

LARGE-SCALE INSTRUMENT MODELS

Large-scale replicas of basic scalers and curettes. Made of durable Immunity Steel®, and mounted on removable display bases. Used to help familiarize students with critical instrument surfaces, designs, and develop a better understanding of their application and maintenance. Models are 12"/30.5cm high.

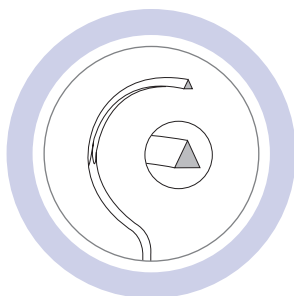

*All large scale models
are also available in
cross section*



#H6
Sickle Scaler

LSMH6

#13
Columbia
Curette

LSM13


#14
Gracey
Curette

LSM14

#30
Jacquette

LSM30

LARGE-SCALE SHARPENING STONE PLASTIC MODELS

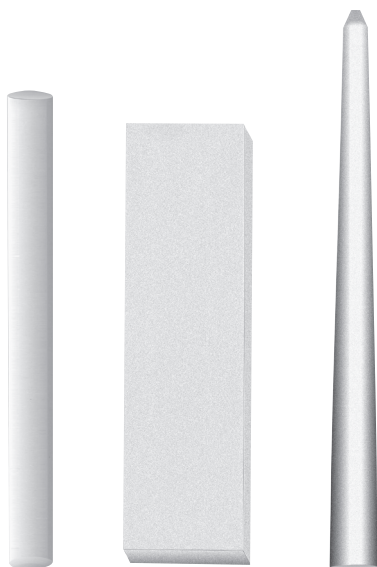

Plastic
Test Stick

LPTS

11.5"/ 29.2cm

#4
Sharpening
Stone


LSS4

3"x 11"x7.5"
(7.6x28x19cm)

#299
Sharpening
Stone

LSS299

11.5"/ 29.2cm

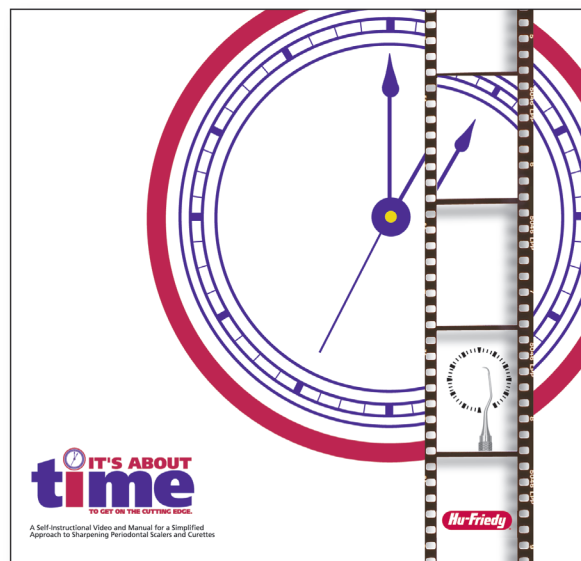

Sharpening
Manual

SHM

It's About Time to Get on the Cutting Edge
Instrument Sharpening manual. Includes post test
for C.E.U. credits.

Sharpening
Video

SDVD


It's About Time to Get on the Cutting Edge
Instrument Sharpening DVD. Manual included with
post test for C.E.U. credits. Self instructional. Time:
12 minutes.



MISCELLANEOUS

SIGNATURE SERIES® GRIPS

Signature Series grips can only be used with Satin Steel Colours®, EverEdge® technology and ResinEight™ instruments, and are used for identifying and organizing instruments. Useful for maintaining proper instrument order in cassettes and identifying instruments by procedure. Will withstand all approved heat sterilization methods.

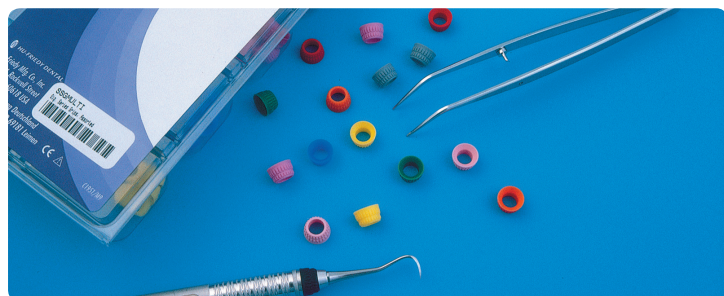

Gray	SSG1	Blue	SSG8
White	SSG2	Black	SSG9
Orange	SSG3	Pink	SSG10
Brown	SSG4	Purple	SSG11
Yellow	SSG5	Ocean Blue	SSG12
Red	SSG6	16 Assorted	SSG0
Green	SSG7	80 Assorted	SSGMULTI

All single color packages contain 16 pieces.
The assorted color packages contain your choice of either 16 or 80 pieces.



Brush Handle
1

HCB1

Used for pit and fissure sealant, acid etch, and placement procedures. Made of Immunity Steel® to allow sterilization by autoclave, chemiclave, and dry heat. Design holds most manufacturers' brush inserts.

Immunity Steel®
Pumice Cup,
Modified

PCUPM



Used to hold prophylaxis paste or medicaments.
Made of corrosion resistant, acid-proof stainless steel.



GF-1
Left

PMGF1

GF-2
Right

PMGF2

Goldman-Fox Periodontal Pocket Marker Right and left paired calibrated (1-2-3-5-7-8-9) designs. Used to establish exterior puncture mark on gingiva at base of pocket to indicate initial line of incision.



Orban Loop

ML

Medicament Loop For subgingival application of medicaments.