

Dome

Endo motor cum Apexlocator

• Build-in ApexLocator

With Multi-frequency length measuring technology the function of root-canal length measuring and enlargement

• 9 Custom Programs


The device can save 9 programs,each program can set different parameters according to the user

• 360° Rotary Angle

The reciprocating motion angle is adjustable FWD-30°60°90°
REV-150°180°210°250°370°

• Colorful OLED Screen

Color OLED displays are available in graphic formats and offer better contrast and viewing angles compared to standard LCDs

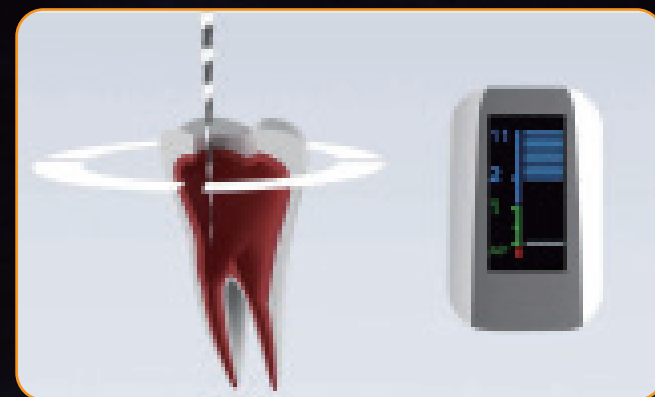

Specifications:

Battery : 3.7V/1500mAh

Rated input voltage : 110v/240v 50/60Hz

output voltage : DC 5V

Revolving speed : 150-1000rpm

Torque : 0.6-4.0N·cm

Angels of FWD-Clockwise : 30°60°90°

Angels of REW-Counter Clockwise :

150°180°210°250°370°